



PTI Wireless Scanner



Less is More





Upgraded chips

The chip upgrade has improved the scanning efficiency of smoke and blood stains, bringing higher quality scanning data.

Detachable tips

You can choose suitable tips according to the different scanning scenarios.



(h)

Start/Pause scanning

Double Click: Switch buccal/Switch color

Long Press: (h) Next step/Model processing

Metal button Remote control of scanning process.

Light weight Only 237g

Built-in battery

Set and charge. No need to replace the battery, making 24/7 standby accessible.

Wireless connection

Breaking free from the constraints of traditional cables, making scanning more convenient and efficient.

3 in 1 The dock The charger The base



No more Tedious Steps 2 H continuous scanning High speed and high-quality transmission



Upper jaw+ Low jaw+ Occlusion Total scan time of 2 min

Super charge 1 H

It only takes 1 hour to reach the full power

Scan up to 60 cases with full battery

Built-in battery

When you are not scanning, you can set it in the base for charging, which greatly extends the battery life



Platform

We have self-designed software to assist you in completing digitalization
1. PTI Center One-stop case management and data scanning
2. PTI OrthoSim Orthodontic simulation contributes to a better doctor-patient communication 3.
Cloud Data storage and fast data transmission




PTI Center

The new software PTICenter can achieve full process connectivity and information integration from case creation, scanning, uploading data, to diagnostic analysis, allowing dentists to enjoy a more convenient and intelligent experience.

	2256	ORDER				1	3.000	20576a	SOFTWARE MA	NAGEMENT	(1) Non service service	6656	ORDE	They be a construction of the			
P = 10 Res Res<	a		**** ****					S. 1997	Installed Software			Ree.					Q01.
1 1 </th <th></th> <th></th> <th></th> <th> 199</th> <th></th> <th>-</th> <th></th> <th>S</th>				 199		-											S
1 State Large 1 State Larg						ACCESS 1			0.8	6 0 x					701944		8:
				-	ALC: U.S. MARCH		2 10 0							ter (sor) and		• (3)	0.
							8 🔳 0										. () ·
			-								du and he						.g.*
APT																STREE STREET	
		-									all he can all		-	(a careche)		10.424	19.02

PTI OrthoSim

The new version of orthodontic simulation software combines advanced technological means and intelligent algorithms to achieve functions such as three-point alignment, intelligent segmentation, and intelligent tooth arrangement, making orthodontic treatment more convenient and accurate, bringing a new experience to both the dentists and the patients.

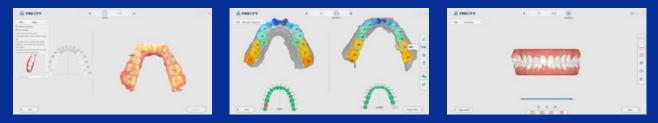

Cloud

PTI scanner provides stable, reliable, and efficient data transmission services for clinics and labs, while also possessing functional features such as case management, 3D preview, and case sharing, helping users better manage and process business data, and improving their work efficiency.





Software Features



Undercut Check

Auxiliary confirmation in place, colour gradient indicates undercut, accurate prediction, and improves the communication efficiency between doctors and technicians.



Margin Line Drawing

Draw accurate and clear margin lines, linking doctors and technicians to reduce rework.



Pre-preparation Scan

Obtain pre-preparation tooth state, retain bite data, and reshape the original good teeth.



Color Switch

True color jaws, clearly restore details, one click to switch colors, and check model data comprehensively.



Al Intelligent Scan

Easily identify soft and hard tissues, eliminate redundant data, and make accurate dental models.



Orthodontic Simulation

Generate simulation animations to help doctors and patients communicate, share quickly, and promote orders.






Bite Check

Color gradient indicates distance, accurately judges the restoration distance and assists in improving the effect of tooth preparation.

Data Lock

Locks the tooth position to prevent data loss and reject excess data interference.

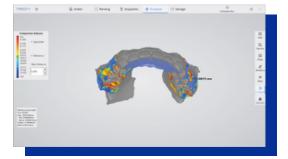

Disease Course Monitoring

Full-cycle diagnosis and treatment monitoring, comparative observation of model differences, simultaneous analysis, comprehensive upgrade of diagnosis and research.



Cuff Tool

Clear cuff data, avoid tissue interference and lock automatically, significantly improving data reliability.

.00

Specification

PRODUCT INFORMATION

Category	PTI Wireless Scanner
Technology	Dynamic Staring 3D scanning
Handpiece Dimension	238mm×42mm×35mm
Handpiece Weight	237g
Tips	6 Tips
Tips Dimension	83x19x14(N) / 82x20x16(D) /
Scan Field	88x20x16(M) 18mm×16mm
Scan Depth	0~15mm
Output Files	STL, PLY, PTY
Connection Method	Wireless
Full Arch Scanning Time	< 2 mins
Light Source	White LED+ Laser
Accuracy	≤15 µm
Precision	≤10 µm
Multi-function Button	YES
Montion Sensing	YES

SOFTWARE

Category	PTI Wireless Scanner
Types of Scanning	Orthodontics, Implant, Restoration
Image Capture	Yes
3D Data Preview	Yes
Output Files	STL, PLY, PTY (Encrypted for user)
Cloud Service	Yes
Cloud Service Price	5GB/Free
Direct Compatibility with Dental CAD	Exocad*
Main Image Processing Tools	Trim/Edit/Recover/Lock, etc.
Margin Line Drawing	Yes
Occlusal Match	Yes
Undercut	Yes
Data Comparison	Yes
AI Technology	Auto Tissue Remove, Overlap Correct, etc.
Software License	Free
Software Upgrade	Free
Transfer Data with API	Available

*All names, logos, brands and other trademarks featured or referred to within this brochure are the property of their respective trademark holders.

PC REQUIREMENTS

Category	PTI Wireless Scanner
CPU	Intel i7-12700H/Intel i7-11800H or
RAM	above 16G/32G
Hard disk	SSD 512G or above
GPU	RTX2060/RTX3060 or above
Operating system	Windows10/11 64bit
Display	Resolution 1920×1080 and above

*Check the PC Recommendation Catalog on our website for more detail.

WARRANTY

Category	PTI Wireless Scanner
Standard Warranty	2 years
Extended Warranty	2 years

*Note: Free maintenance during the first two years (No labor fee charged)

Customer Support



www.ptiscanner.com







PTI SCANNER ALL OVER THE WORLD